


The exclusive advantages of the products



Comfort Zone

Soft zones in the sensitive flexor areas ensure that the support does not press on or cut into the skin, even when worn all day long.

> The result: soothing pain relief that simply feels good.

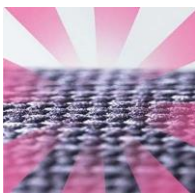

Clima Fresh

Clima Fresh has antibacterial properties built into the knit of the fabric. This prevents bacteria and micro organism from colonising in the fabric or replicating, thus preventing unpleasant odours.

Clima Fresh



Clima Comfort

The elaborately worked knit ensures a balanced microclimate of the skin. Moisture is transported reliably and quickly from the inner surface to evaporate from the outer surface.

At the same time, the material's active breathing structure guarantees the necessary oxygen supply.

> The result: patients don't work up a sweat with their support and find it comfortable – even in summer.



Airtex+

The elastic spacer fabric is the combination of medical compression and the textile effects **SMART TEMP** and **FRESH TECH**.

SMART TEMP

The adaptive thermoregulation is based on an innovative polymer that changes its properties in response to temperature. At low temperature the polymer binds moisture in its structure keeping the body drier. At high temperature it releases the stored moisture and cools the textile and skin.

FRESH TECH

The antibacterial effect of FRESH TECH guarantees deep optimal hygienic conditions even after long term use.

> The combination of medical compression and functional textile effects guarantees improved wearing comfort.



Airtex Material

An innovative, breathable, light and compressive material for comfortable activation of proprioception. It optimises moisture transport and ensures the highest comfort, even in hot climates.

Spinal Orthosis

Lumbamed® disc

Lumbar support
with adjustable
stabilisation system



Indications

- intervertebral disc protrusion or disc prolaps (nonsurgical, postoperative)
- recurring lumbar dysfunction
- spondyloarthrosis
- lumbar overload syndrome
- static muscular imbalance (lumbar and thoracolumbar area)
- all indications requiring reduction of lordosis



Mode of action

- stabilisation of the lumbar spine area and the lumbosacral transition by means of integrated dorsal stabilisation elements
- reduction of the lumbar lordosis via abdominal compression, strap system additionally supports the abdominal area
- due to the modular design (back pad, straps and exchangeable stabilising rods) the stability of the support can be adapted / reduced to suit the course of recovery

Product benefits



Diamond-shaped functional zone
Large air holes in the back allow for increased breathability.

silver



Clima Comfort



The support is particularly simple to apply thanks to cleverly positioned hand loops and improved elasticity of the comfort knit.



Clima Fresh




Stability can be adapted in three stages to suit the treatment (pad with bridging effect incl. removable tension straps, interchangeable hardness of the back rods).




Comfort Zone

women

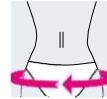

circumference hip (cm)	76 - 86	86 - 98	98 - 111	111 - 125	125 - 140	140 - 156
size	I	II	III	IV	V	VI
article	K.220.3	K.220.3	K.220.3	K.220.3	K.220.3	K.220.34

men

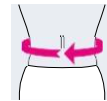

circumference waist (cm)	66 - 76	76 - 88	88 - 101	101 - 115	115 - 130	130 - 146
size	I	II	III	IV	V	VI
article	K.220.33	K.220.33	K.220.33	K.220.33	K.220.33	K.220.33

Lumbamed® basic

Lumbar support with
Clima Comfort fabric



Indications silver

- lower back pain (lumbago), ligamentous disorders
- myostatic, muscular insufficiency
- irritative conditions of the facet joints
- symphysiolysis
- postdissectomy syndrome
- irritative conditions of the lumbosacral transition

Mode of action

- stabilisation of the lumbar spine area and the lumbosacral transition by means of 4 integrated dorsal spring ribs
- reduction of lordosis of the lumbar spine by means of proprioceptive compression in the abdominal area

Product benefits



The support can be opened comfortably with a single finger (due to the triangular tab).

silver



Clima Comfort



The fastener can be adapted continually to the abdominal contours.

The edges "bend out" comfortably - especially when sitting (stops them pressing into and chafing the skin).



Clima Fresh



The support is particularly simple to apply thanks to cleverly positioned hand loops and improved elasticity of the comfort knit.




Comfort Zone

women

men

I Orthosis

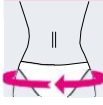

circumference hip (cm)	76 - 86	86 - 98	98 - 111	111 - 125	125 - 140
size	I	II	III	IV	V
article number	K.200.Δ4	K.200.Δ4	K.200.Δ4	K.200.Δ4	K.200.Δ4



circumference hip (cm)	66 - 76	76 - 88	88 - 101	101 - 115	115 - 130
size	I	II	III	IV	V
article number	K.200.Δ3	K.200.Δ3	K.200.Δ3	K.200.Δ3	K.200.Δ3

Shoulder Orthosis

medi SAS® comfort

Shoulder abduction
support; 15° abduction



Indications

- after rotator cuff repair
- after shoulder stabilization
- after first-time anterior shoulder dislocation (nonsurgical)
- fractures of the humerus head (nonsurgical, postoperative)
- after insertion of shoulder prosthesis
- postoperative treatment after acromioclavicular jointluxation
- biceps tendon tenodesis

Mode of action

- lateral positioning in 15° abduction with anatomicallyshaped cushion to avoid
- relieves pressure on the subacromial space
- supports ventral limb and capsule areas

Product benefits



Handrest (with integrated exercise ball) enables patients to exercise the forearm and wrist

- improved circulation
- prevents stiffening of the wrist.




Anatomically shaped cushion guarantees comfortable arm support.



The special ring-strap-solution transmit tension forces take away from the healthy shoulder / neck region

silver

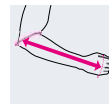

length of forearm (cm)	up to 34	34 - 39	over 39
size	S	M	L
article number	R.013.01	R.013.01	R.013.01

◊ = fill in for orientation: 2 = left, 3 = right

Omomed

Shoulder support with
movement-limiting function



Indications

- proximal humeral head, scapula or glenoid fractures andAC joint dislocations (nonsurgical, postoperative)
- nonsurgical with anterior shoulder dislocation and chronic anterior subluxation as well as painful omarthrosis
- postoperative after shoulder stabilisations, rotator interval reconstruction, refixation of the biceps tendon anchor

Mode of action

- limitation and prevention of injurious movements thanks to the adjustable limitation of external rotation, abduction and anteflexion of the shoulder joint
- the elastic knit exerts gentle local compression on the superficial soft-tissues to limit and counteract swelling and reduce muscular tension
- corrects postural deficits in the shoulder with accurate centring of the head of the humerus in the glenoid cavity and easing of shoulder pain

Product benefits



Flexible strap system permits individual adjustment and variable limitation of movement at the shoulder joint.



Intelligent strap positioning ensures comfortable wearing and unrestricted function, even in everyday clothing.

silver



Clima Comfort



Clima Fresh



Comfort Zone



circumference upper arm (cm)	up to 26	26 - 29	29 - 33	33 - 36	36 - 39
size	0	I	II	III	IV
article	K.300.20	K.300.201	K.300.20	K.300.203	K.300.20

◊ = fill in for orientation: 2 = left, 3 = right

Shoulder Orthosis

medi Epico ROMs

Elbow brace for mobilisation with immobilisation of the PRU joint



Indications

- treatment of elbow dislocation with associated injuries (nonsurgical / postoperative)
- treatment of medial or lateral epicondyle fractures
- treatment of fractures of the radial head (nonsurgical / postoperative)
- treatment of distal humeral fractures (nonsurgical / postoperative)
- postoperative care after elbow arthroplasty

Mode of action

- immobilisation of the elbow joint by fixation of the upper and lower arm
- joint movement is limited to 0° - 120° to guarantee ROM appropriate to the stage of rehabilitation
- prevention of movements that can cause luxation or subluxation
- prevention of pronation and supination movement due to hand rest

Product benefits



Limitation 0° - 120° serves to immobilise the elbow joint by fixation of areas of the upper and lower arm.



Flat Quick Set hinge pleasant to wear, even under clothing. The hinge can be simply adjusted by pressing a button.



Padding made of innovative foam and fits tightly and reliably – the fluted surface slightly compresses the muscle receptors, stimulating metabolism and preventing muscle tension.



Aluminium construction can be precisely shaped and adapted to the shape of any arm.

Length of the hand cushion can be adjusted – the orthosis is securely attached to the arm.



length of forearm (cm)	up to 29	29 - 42
size	S	L
article number	R.032.00	R.032.004
◊ = fill in for orientation: 2 = left, 3 = right		

grey / graphit

Hip Orthosis

medi Hip one

Functional hip joint orthosis



Product benefits



Hip parts are easy to bend.



Comfort due to cushioned padding.

black / white

- dynamic abduction

Indications

- nonsurgical hip therapy and dislocation of the hip
- traumatic dislocation of the hip joint
- acetabular fracture
- postoperative conditions
- hip prosthesis displacement
- replacement of hip prosthesis
- acetabular reconstruction
- resection of the head of the femur
- adjustment osteotomy



Mode of action

- guards against dislocation / luxation
- controlled early functional exercise therapy
- extension and flexion limitations set in 15° increments up to 90°

Variable angle adjustment.



(A) circumference hip part (cm)	70 - 90	90 - 110	110 - 135
size	S	M	L
article number	5.5F9.102	5.5F9.103	5.5F9.104


◊ = fill in for orientation: 2 = left, 3 = right

(B) circumference thigh part (cm)	38 - 50	50 - 62	62 - 76
size	S	M	L
article number	5.5F9.202	5.5F9.203	5.5F9.104

◊ = fill in for orientation: 2 = left, 3 = right

To get the complete product please use two article numbers: hip part and thigh part measurement : A circumference hip

B circumference thigh measured 5 cm below the crotch

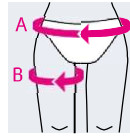

Knee Orthosis - Soft Ligament Braces

Genumedi®

Knee support with patella silicone ring



Indications

- slight instability of the knee joint
- articular effusions and swellings
- chronic, postoperative and post-traumatic inflammation
- osteoarthritis and arthritis
- patello femoral pain syndrome

Mode of action

- massage effect silicone insert
- improves the blood circulation and aids absorption of haematomas and oedemas
- improves proprioception and muscular support (stabilisation of the knee)



Product benefits



Integrated silicone patella donut with 3D profile for a gentle massaging effect.



Special inner border for the silicone ring provides maximum pressure relief for the patella.



Integrated silicone dots ensure an effective grip

silver



Clima Comfort




Clima Fresh



Comfort Zone

Genumedi®

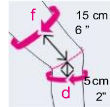

circumference "f" (cm)	37 - 40	40 - 43	43 - 46	46 - 49	49 - 52	52 - 55	55 - 58
circumference "d" (cm)	25 - 28	28 - 31	31 - 34	34 - 37	37 - 40	40 - 43	43 - 46
size	I	II	III	IV	V	VI	VII

Genumedi® pro

Orthosis for physiological guidance of the knee joint



Indications

- mild to moderate instability of the knee joint
- osteoarthritis of the knee
- mild meniscus injuries
- mild collateral ligament instabilities
- postoperative stabilisation and compression, for example following partial meniscectomy or cartilage smoothing

Mode of action

- massage effect from the integrated silicone pad
 - reduces swelling (oedemas, haematomas)
- proprioceptive compression
 - guides and stabilises the joints
- easyglide hinges closely mimic the natural roll and glide motion of the knee joint
 - providing for comfortable stabilisation

Product benefits



Modern easyglide hinges stabilise and support the natural roll and glide motion of the knee joint.



With additional straps for stabilisation and ultimate wearing comfort.



Integrated silicone patella donut with 3D profile for a gentle massaging effect.



Integrated silicone dots ensure an effective grip

silver



Clima Comfort




Clima Fresh



Comfort Zone

Genumedi® pro

Knee Orthosis - Soft Ligament Braces

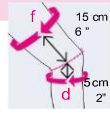

circumference "f" (cm)	37 - 40	40 - 43	43 - 46	46 - 49	49 - 52	52 - 55	55 - 58
circumference "d" (cm)	25 - 28	28 - 31	31 - 34	34 - 37	37 - 40	40 - 43	43 - 46
size	I	II	III	IV	V	VI	VII

Knee Orthosis - Soft Ligament Braces

Stabimed®

Short, functional soft brace with flexion / extension limitation



Indications

- mild to moderate instability, chronic instability
- meniscus injuries
- traumatic or degenerative changes at the knee joint
- mild to moderate osteoarthritis

Mode of action

- safe collateral stabilisation of the knee joint
- extension and flexion limitation protects against damaging joint movements and relieves injured structures
- improvement of proprioception to increase the own guidance and stabilisation of knee joint
- prevention of hyperextension
 - extension: 0°, 10°, 20°, 30°
 - flexion: 0°, 10°, 20°, 30°, 45°, 60°, 75°, 90°
- immobilisation: 0°, 10°, 20°, 30°

Product benefits



Very comfortable wearing due to persistent net material on the back.


grey



Tool-free wedge replacement.



medi Airtex

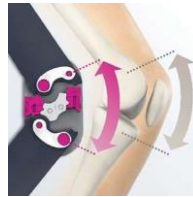

The physioglide hinge mimics the rolling / gliding movement of the knee and ensures a secure fit.



wrap version

circumference thigh (cm)	31 - 37	37 - 45	45 - 51	51 - 57	57 - 65	65 - 74	74 - 84
size	XS	S	M	L	XL	XXL	XXX L



Following wedges are pre-installed on the brace at the point of delivery: extension 0° / flexion 90°.

medi PT control

Soft brace for patella tracking control with flexion / extension limitation



Indications

- nonsurgical or postoperative treatment after patella luxation with mild knee instability
- lateral malalignment with mild knee instability
- lateral release surgery
- retropatellar cartilage surgery

Mode of action

- reliable relief of injured or operative treated structures
- correction of lateral Patella maltracking
- extension and flexion limitation protects against damaging joint movements and relieves injured structures
- limitations:
 - extension limitation at: 0°, 10°, 20°, 30°, 45°
 - flexion limitation at: 0°, 10°, 20°, 30°, 45°, 60°, 75°, 90°
 - immobilisation at: 0°, 10°, 20°, 30°, 45°
- collateral stabilisation of the knee joint

Product benefits



Tool-free wedge replacement

The wedges can be changed quickly and easily with built-in tweezers with no need for screws.

grey



Also available as wrap version

Excellent compliance due to wrap-around version with hook and loop fastener closing to the front.




medi Airtex



Correction of patella glide due to an individually adjustable patella stabilisation buttress and air-cushion.




wrap version

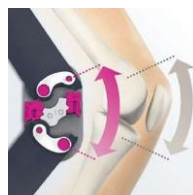

The physioglide hinge mimics the rolling / gliding movement of the knee and ensures a secure fit.

Knee Orthosis - Soft Ligament Braces

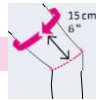

circumference thigh (cm)	31 - 37	37 - 45	45 - 51	51 - 57	57 - 65	65 - 74
size	XS	S	M	L	XL	XXL

□ = fill in variant: 2= medi PT control, 3= medi PT control wrap / ∅ = fill in orientation: 2= left, 3= right
 Following wedges are pre-installed on the brace at the point of delivery, extension 0° / flexion free

M.4s[®] PCL dynamic

4-point rigid frame brace with extension/ flexion limitation for treating injuries of the posterior cruciate ligament



Indications black

- nonsurgical treatment after injury of the posterior cruciate ligament
- postoperative treatment after reconstruction of the posterior cruciate ligament

Mode of action

- the protection pad actively exerts a force directed anteriorly to hold the tibia in the desired position
- individual and accurate adjustment of the protection pad due to integrated tensioning dial
- safe collateral stabilisation by rigid frame
- extension and flexion limitation protects against damaging joint movements and relieves injured structures
 - extension limitation at: 0°, 10°, 20°, 30°, 45°
 - flexion limitation at: 0°, 10°, 20°, 30°, 45°, 60°, 75°, 90°
 - immobilisation at: 0°, 10°, 20°, 30°, 45°

Product benefits


black



PCL protection pad
the pad presses actively on the tibia and holds it in the desired position.




Tensioning Dial
the integrated tensioning dial enables the protection pad to be adjusted exactly to suit the individual patient.

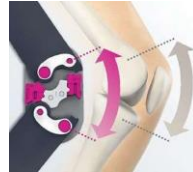

The physioglide hinge
mimics the rolling / gliding movement of the knee and ensures a secure fit.



Pre-shaped tibia pad
the winged shape of the tibia cushion moulds itself perfectly to the anatomical shape of the shin bone.

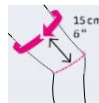

circumference thigh (cm)	31 - 37	37 - 45	45 - 51	51 - 57	57 - 65	65 - 74
size	XS	S	M	L	XL	XXL

article number G.08100 G.08100 G.08100 G.08100 G.08100 G.08100
 ∇ = fill in for orientation: 2 = left, 3 = right
 Following wedges are pre-installed on the brace at the point of delivery: extension 0° / flexion 90°



M.4s® OA comfort

4-point knee orthosis for varus or valgus pressure relief and stabilisation




Indications

- severe medial or lateral gonarthrosis (at least level III) combined with ligament instability
- unicompartmental relief before HTO (brace test)
- for relief of painful stress fractures (e.g. tibiahead)
- cartilage surgery / chondroplasty / meniscus refixation / displacement osteotomy with required postoperative relief of the medial or lateral compartment
- strong overloading of the medial or lateral compartment
- permanent relief, where surgery is contraindicated

Mode of action

- relief of the lateral or medial compartment using the 3-points principle
- pain relief by reducing weight-bearing on the affected side of the joint
- optimal stabilisation using 4-point principle
- safe collateral stabilisation by rigid Frame
- extension and flexion limitation protects against damaging joint movements and relieves injured structures
- limitations:
 - extension: 0°, 10°, 20°, 30°, 45°
 - flexion: 0°, 10°, 20°, 30°, 45°, 60°, 75°, 90°
 - immobilisation: 0°, 10°, 20°, 30°, 45°

Product benefits

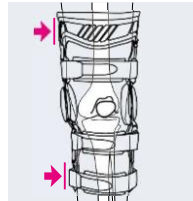

3-point principle for guaranteed protection.



Movable thigh clamp for a perfect contact and positioning on the leg.




Double-layer condylar padding cover to increase comfort and reduce pressure friction.



Ultra-flat and lightweight frame and hinge construction.

Anatomically preformed and torsion-resistant frame geometry

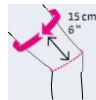

circumference thigh (cm)	31 - 37	37 - 45	45 - 51	51 - 57	57 - 65	65 - 74
size	XS	S	M	L	XL	XXL
article number	G.52.0.001	G.52.0.002	G.52.0.003	G.52.0.004	G.52.0.005	G.52.0.006

□ = fill in for variant: 7 = Varus, 8 = Valgus / ◊ = fill in for orientation: 2 = left, 3 = right



Knee Orthosis - Post-Op Braces

medi ROM / medi ROM cool

Universal knee splint with flexion /
extension limitation and options for
immobilisation



Product benefits



Patented "Quick Set"
hinge.

Easy, comfortable 'click
closing' strap fasteners.

Adjustable skin friendly
padding.

grey

Indications

- postoperative immobilisation and for later mobilisation
- for secure stabilisation of injured or operated collateral ligaments

Mode of action

- collateral stabilisation of the knee joint
- postoperative immobilisation and controlled mobilisation due to variable extension and flexion limitation controlled mobilisation
- limitation in steps of 10° to 15°
 - flexion: 0°, 10°, 20°, 30°, 45°, 60°, 75°, 90°, 105°, 120°
 - extension: 0°, 10°, 20°, 30°
 - immobilisation: 0°, 10°, 20°, 30°

medi Classic

Knee immobilisation brace at
0° or 20°



Product benefits

Individual, rapid adjustment to fit varying leg
circumferences.

Removable medial, lateral and posterior
stabilisation rods.

Comfortable, washable, skin friendly terry
towelings cloth material.

Available as 0° or 20° version

blue



medi Kidz –
also available as
Kidz-Version.

Indications

- postoperative immobilisation
- preoperative pain conditions
- conditions following patella dislocations

Mode of action

- immobilisation of the knee at 0°
resp. 20°



Ankle Orthosis

Levamed

Ankle support with silicone cushions



Indications

- slight instability of the ankle joint
- sprains and distortions
- joint effusions and swellings due to arthritis and osteoarthritis
- chronic, postoperative and posttraumatic inflammation

Mode of action

- massage effect due to integrated silicone pads
→ reduction of swellings (oedema, haematomas)
- proprioceptive compression
→ guidance and stabilisation of the joints



circumference ankle (cm)	17 - 20	20 - 22	22 - 24	24 - 26	26 - 28	28 - 30
size	I	II	III	IV	V	VI
article number	K.010.Δ1 1	K.010.Δ1 2	K.010.Δ1 3	K.010.Δ1 4	K.010.Δ1 5	K.010.Δ1 6

Δ= fill in for colour: 2 = silver

Product benefits



Innovative comfort knit with anatomically tailored support, soft instep and supple edging.

silver



Clima Comfort



Integrated silicone pad with a 3D profile for a gentle massage effect and to reduce swelling.



Clima Fresh

Accessories



A donning aid / slip-in sock is supplied with each product.



Comfort Zone

Ankle Orthosis



Ankle Orthosis

Levamed active

Ankle support with stabilisation strap system



Indications

- joint sprain with overstretch of the capsule
- chronic ligamentous instability with recurrent supination traumas (twisting of the foot)
- in the second phase of treatment for deltoid ligament injuries
- acute or chronic post-traumatic or postoperative irritation of the talocrural and subtalar joints
- nonsurgical treatment of injuries of the syndesmosis not requiring surgery
- follow-up treatment of fractures not requiring a cast

Mode of action

- individually adjustable and removable strap system with elastic and inelastic sections
→ additional hold and stabilisation, protection "against twisting (supination)"
- massage effect by the integrated silicone pads
→ reduction of swelling (oedema, haematomas)
- proprioceptive compression
→ guidance and stabilisation of the joints

Product benefits



Clever belt design and intuitive belt guidance for easier donning.



Individually adjustable and removable belt system with elastic and non-elastic areas.




Easier donning and doffing thanks to the V-stretching zone at the upper cuff.

Accessories



Accessories Easy donning with the enclosed balloon silk donning aid.

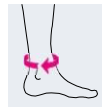


circumference ankle (cm)	17 - 20	20 - 22	22 - 24	24 - 26	26 - 28	28 - 30
size	I	II	III	IV	V	VI
article number	K.011.2	K.011.2	K.011.2	K.011.2	K.011.2	K.011.2
◊ = fill in for orientation : 2 = left, 3 = right						

silver



Clima Comfort

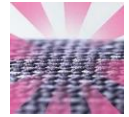

Clima Fresh



Comfort Zone

Levamed stabili-tri

Triple-action modular orthosis for the functional care of ankle injuries



Product benefits



Especially wide opening tongue for comfortable donning and doffing.

rapid absorption of haematomas and oedemas

blue / azure



Clima Comfort



Indications

- ankle sprains (grades II to III)
- chronic instability of the talocrural and subtalar joints
- postoperative care after surgery on the ligamentous apparatus of the ankle joint


Mode of action

- modular design supports functional treatment
- adaptable to the three healing phases - for safe and reliable treatment and mobilisation of the patient
- talus lock tongue reduces anterior talar shift
- medial stabilising element gives support and relieves stresses on the lateral side of the foot
- stabilising strap guides the foot in gentle pronation and supports the healing process
- compressive knit improves proprioception and promotes more

Ankle Orthosis

oss- lacing offers flexible pressure customisation and simple handling thanks to click function-laces 'click' quickly and reliably into the eyes.

The medial position of the stabilising element protects the injured lateral side of the foot and ensures high comfort in wear.

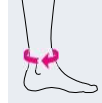

circumference ankle (cm)	17 - 21	21 - 24	24 - 27	27 - 30
size	I	II	III	IV
article number	K.012.30	K.012.30	K.012.30	K.012.40

◊ = fill in for orientation : 2 = left, 3 = right

Ankle Orthosis

M.step®

Functional ankle brace
with foam gel padding



Product benefits

 black



Only two straps allow a quick and secure fastening; straps can be taken off the shells and individually fixed.

Very light foam gel padding with soft, skin friendly surface.

Indications

- nonsurgical and postoperative treatment of fibular ligament ankle injuries (grades II to III)
- after sprains
- chronic instabilities of the ankle



Mode of action

- rigid shells provide very good stabilisation in pronation and supination movements
- the foam-gel of the padding adapts perfectly to the varying anatomic conditions and provides a uniform compression

medi ROM Walker

Walker with range of motion control



Product benefits

 silver



Patented "QuickSet" hinge.

Indications

- postoperative and nonsurgical treatment of fractures of the fibula
- dislocated fractures of the talocrural joint
- nonsurgical treatment of injuries of the fore foot, middle foot or tarsus
- ligament refixation and peroneal tendon dislocation
- follow-up treatment of fractures of the talus, calcaneus and medial malleolus
- early functional therapy after rupture of the Achilles tendon (nonsurgical / postoperative)

Mode of action

- optional fixed joint angle or partial excursion; this permits early functional mobilisation
- plantarflexion: 0°, 10°, 20°, 30°, 40°
- dorsiflexion: 20°, 10°, 0°
- immobilisation at each 10° angle setting physiological gait due to special sole design
- controlled compression by anterior and posterior shell as well as circular hook and loop fastener straps



shoe size EU	32 - 38	38 - 43	43 - 48
shoe size UK	1 - 5	5 - 9	9 - 13
size	S	M	L
article number	G.160.00 2	G.160.00 3	G.160.00 4

Ankle Orthosis

protect. CAT Walker

Comfortable walker for ankle and foot injuries - in a lightweight shell construction



Indications

- achilles tendon rupture, ligament / soft part and tendon injuries (postoperative)
- metatarsal fractures, achilles tendon injuries, severe distortions (nonsurgical)
- stable foot and joint fractures, fracture postoperative care (talus, luxation of the upper ankle joint)
- ligament refixation

Mode of action

- safe and effective immobilisation
- four strap system ensures high stability and reduces anterior talar shift

Product benefits



Air system with medial and lateral air passage around the ankle and up to the forefoot.



Air pump and ventilation valve for easy handling.



Easy adjustment and removal of toe protector.



Rounded sole for natural, stable gait.

shoe size EU	36 - 40	40 - 44	44 - 47
size	S	M	L
article number	G.950.01 2	G.950.01 3	G.950.014

grey

protect.CAT Walker short

Comfortable walker for ankle and foot injuries - in a lightweight shell construction



Indications

- ligament / soft part and tendon injuries (postoperative)
- metatarsal fractures, severe distortions (nonsurgical)
- stable foot and joint fractures, fracture postoperative care (talus, luxation of the upper ankle joint)
- ligament refixation

Mode of action

- safe and effective immobilisation
- three strap system ensures high stability and reduces anterior talar shift

Product benefits



Air system with medial and lateral air passage around the ankle and up to the forefoot.



Air pump and ventilation valve for easy handling.



Easy adjustment and removal of toe protector.



Rounded sole for natural, stable gait.

shoe size EU	36 - 40	40 - 44	44 - 47
size	S	M	L
article number	G.950.11 2	G.950.11 3	G.950.114

grey